

stryker[®]

Spike and Go with StrykeFlow[™]

StrykeFlow[™] II Suction Irrigation



StrykeFlow™ II Suction Irrigation

Experience the benefits of optimal suction irrigation performance. This completely disposable battery powered suction irrigator includes a small, powerful pump that spikes directly into the bag and achieves exceptional flow rates. Hassle free set-up equates to decreased OR preparation and turnover time. Add ease and simplicity to any laparoscopic procedure with StrykeFlow™ II.

Product Specifications

Maximum Flow Rate of 4 L/min

620 mmHg Pressure

Powered by 8 AA Batteries

5mm Instrument Portal

Features & Benefits


Direct Spike into Bag
for Easy Set-Up

Ergonomic Handpiece
for Ambidextrous Grip and Comfort

Pop-Top Battery Housing
for Easy Battery Removal

Variable Suction Control Valve
for Smoke Evacuation

Compatible with any Stryker Suction
Irrigation Tip

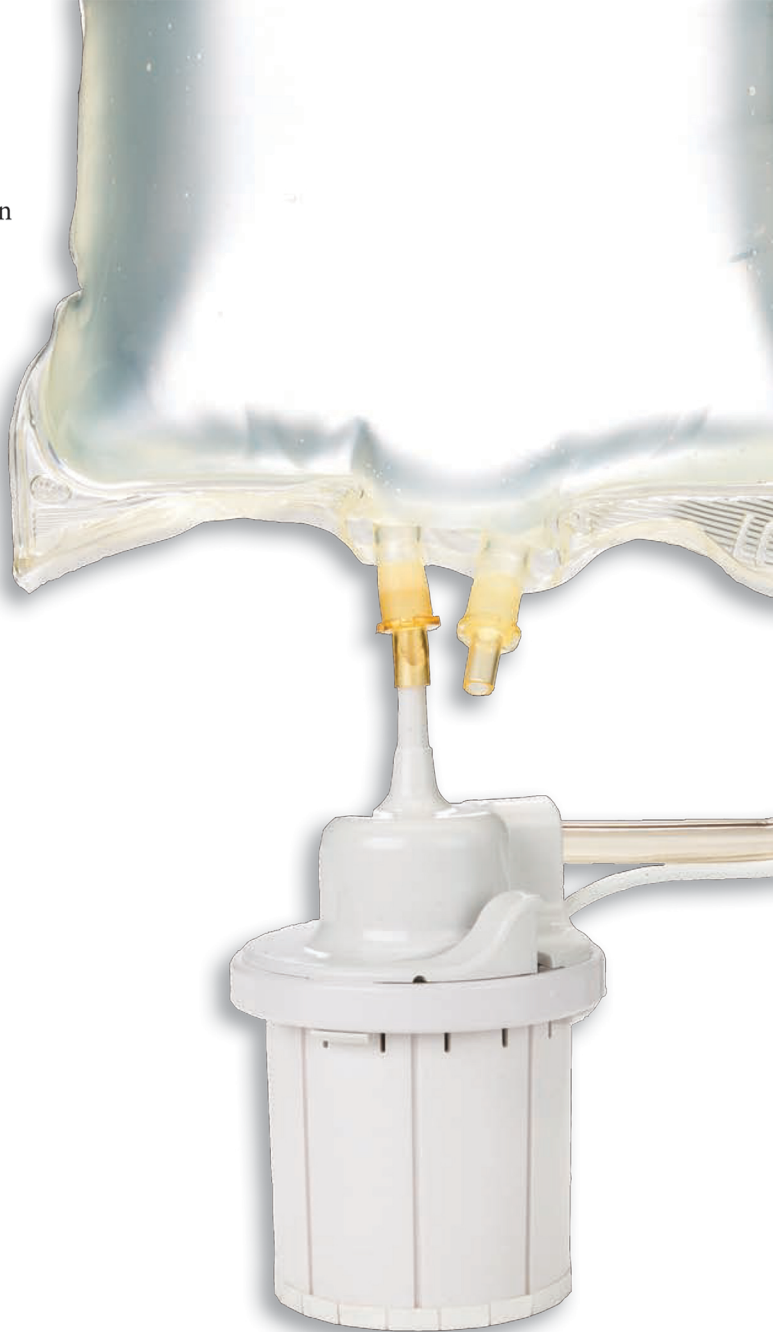

Product Ordering Information

ORDERING INFO	PRODUCT
250-070-500	StrykeFlow™ II No Tip
250-070-520	StrykeFlow™ II with Disposable Tip



EASMED

Tel: (65) 6336 0806
Fax: (65) 6471 3233
customer@easmed.com