



# YIGAO MEDICAL

## UROLOGY PRODUCTS

# ELEPHANT—II

## Disposable Flexible and Negative Suction Ureteral Access Sheath


Maximize Stone-Removal Efficiency

Reduce Intra - Renal Pressure

Clear Surgical Field of View

Avoid Thermal Damage

Reduce Complications

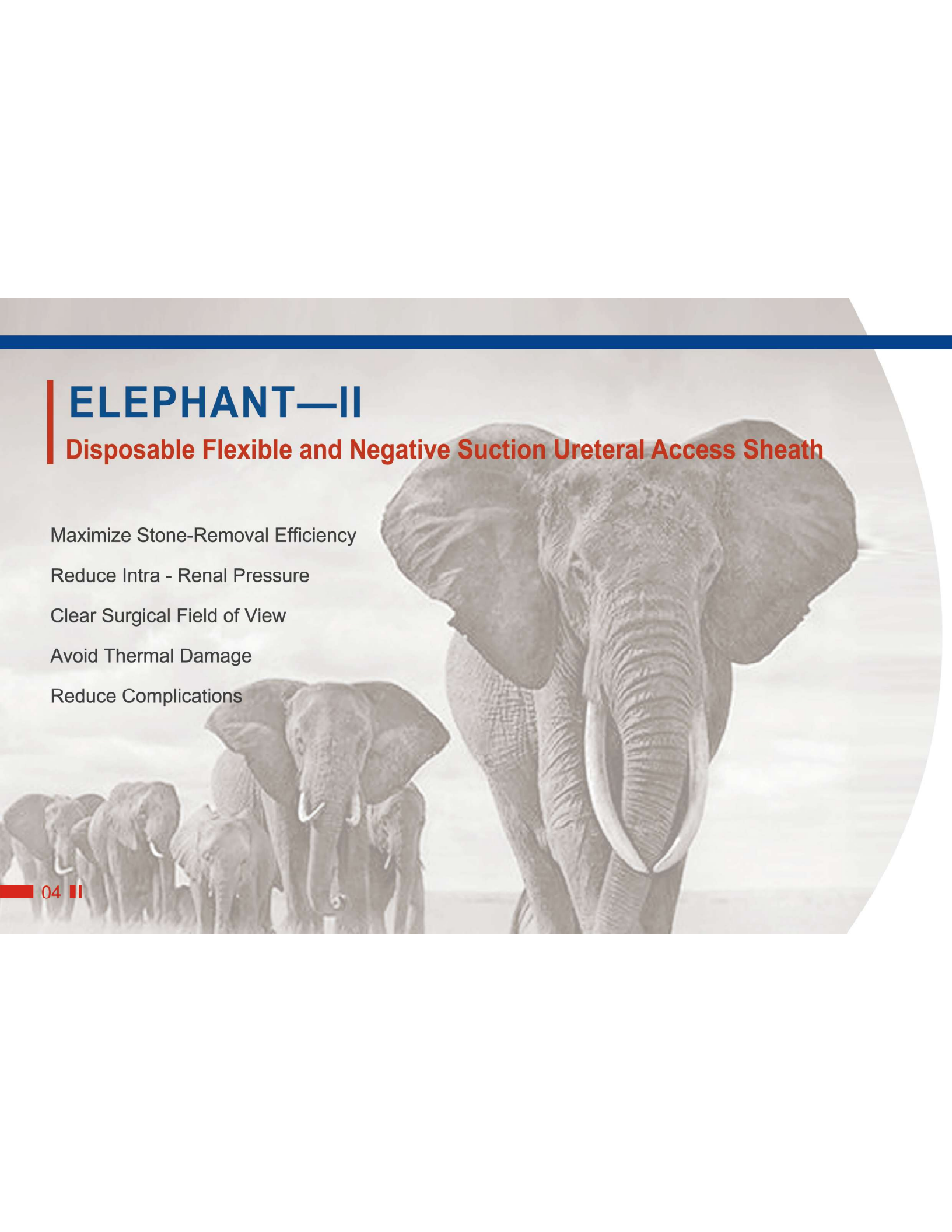

## Dilation

### Specification and Model

Inner Diameter (Fr)	9.5	10	10.5	11	12	12.5	13
Length(cm)	40	50	52	55			

Please contact us for more information (See back cover)




The outer sheath and the 15cm dilator tip are covered with hydrophilic coating for providing ease during insertion.

The new intelligent pressure control module eliminates the need to use the thumb.




Stone Collector  
Used for collecting stones during surgery

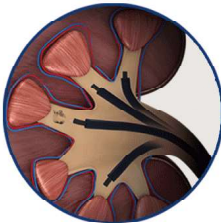

The access sheath can be bent with the ureteroscope to enter the renal calyces and reach the stones to maximize stone - removal efficiency.

## Commence Suction While Stones Are Being Fragmented

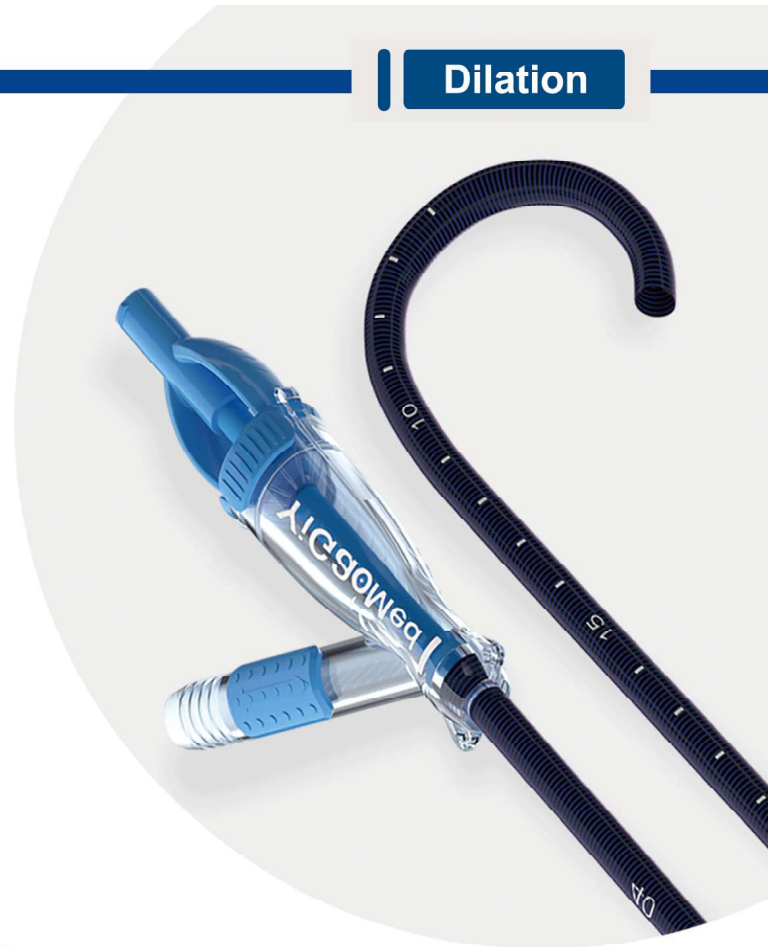

## Dilation

### Specification and Model

Inner diameter(Fr)	9.5	10	10.5	11	11.5	12	12.5	13
Length(cm)	35		45					

Please contact us for more information (See back cover)

The outer sheath and the **15cm** dilator tip are covered with hydrophilic coating for providing ease during insertion.

A spring support ring covers the inner layer of the sheath tube, effectively supporting the working channel and providing an excellent endoscope-to-sheath ratio.

The soft tapered tip of the dilator can make it easier to pass through the twisted parts of the ureter and minimize trauma during insertion.

## Disposable Ureteral Access Sheath

Used to establish a conduit during endoscopic urological procedures, facilitating the passage of endoscopes and other instruments into the urinary tract.

